





## Treatment Result


### • Pigmentation

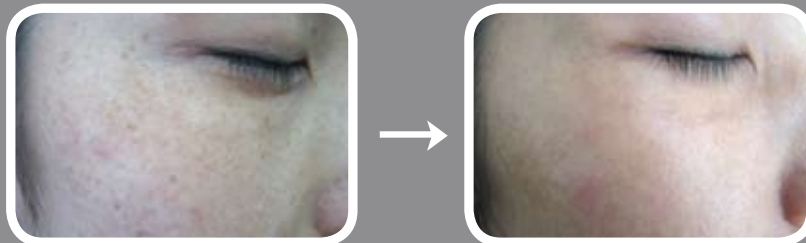

### • Vascular

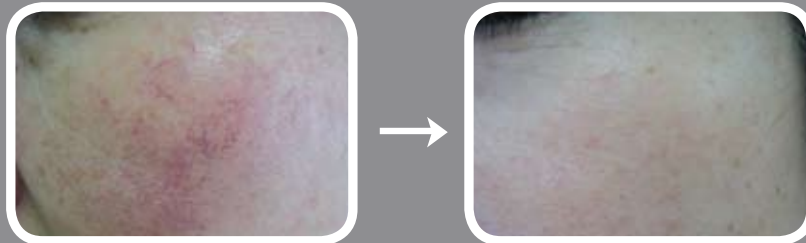

### • Wrinkle

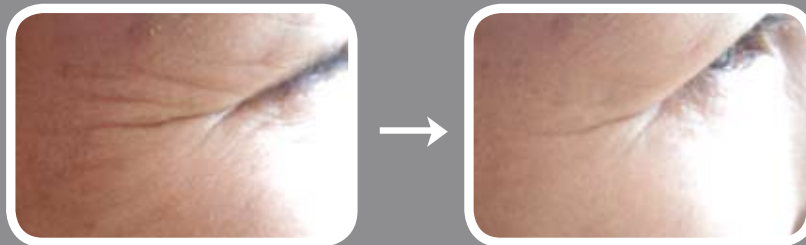

### • Acne



### • Hair removal



Before

After

# Hironic

## Specification

Lamp Type	Xenon Flash Lamp
Wave Length	450-950 nm
Standard Filter	450, 560, 640, SHR Filter
Energy Density	1Pulse 5 ~ 40 J/cm2
Spot Size	10 X 38 mm
Pulse Sequence	Single / Double / Triple pulses
Pulse Duration	5-30 ms
Pulse Delay	5.0-99.0 ms
Repetition Rate	Up to 1Hz
Cooling System	Peltier Cooling System
Control	LG 10.4 Inch LCD Touch Screen
Electrical Requirements	AC220V / 50-60Hz / 30A
Dimensions(W X D X H)	350 X 520 X 1170
Weight	60kg



## Medical Success Partner Hironic



### What's Hironic?

Hironic Co., Ltd. has set a new standard of aesthetic medical devices with our technology and field experiences in medical industry, and pursuing a goal to manufacture the most qualified medical machines. Hironic Co., Ltd. who considers customer`s needs as a top value based on the leading technology and production knowhow would be always with your success as a specialized medical pioneer.



### Hironic Co., Ltd.

#913 Sicox Tower, 513-14, Sangdaewon-Dong, Joongwon-Gu, Sunghnam, 462-120, Korea  
**Tel** +82-31-777-9711(Domestic), +82-70-7633-9582(Overseas)  
**Fax** +82-31-777-9713 [www.hironic.com](http://www.hironic.com)

### Location of the dealer

[www.hironic.com](http://www.hironic.com)



## Safe & Effective High speed & Large tip size

- Hair removal
- Photo Rejuvenation
- Vascular Lesions
- Pigment Lesions
- Acne



## Intelligent Pulsed Light

# MPL





- Hair removal
- Photo Rejuvenation
- Vascular Lesions
- Pigment Lesions
- Acne




Intelligent Pulsed Light

**MIPL**


### 3rd Generation of IPL system

Our concern for patient safety led us to develop this unique Dual Mode Filtering technique, which filters out potentially harmful wavelengths of light. MIPL utilises only relevant wavelengths that are not absorbed by water in the skin tissue. Dual Mode Filtering makes treatments safer, as less energy is needed.

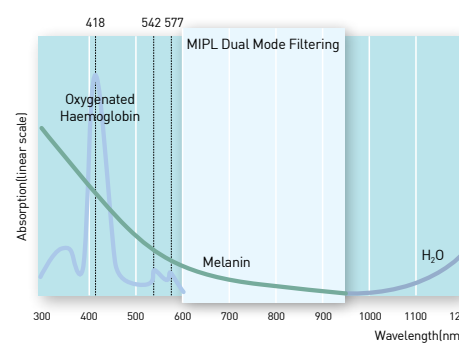

### Square Pulse Technology ensures optimal safety and effectiveness

All IPL systems work by emitting intense light in short pulses. The curved line represents the pulse produced by many IPL systems. But MIPL's Square Pulse Technology ensures that the light goes to the desired level immediately, remains at this level for the duration of the pulse, and ends instantly. As a result, MIPL treatment works optimally and minimises the risk of side effects.

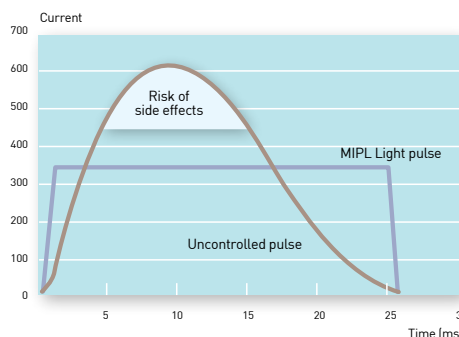

### 1. Specialized Treatment System

Using Dual Mode Filtering technology, patients expect safe and effective treatment.

### 2. Convenient Maintenance and Repair

Separated body and hand-piece, easy to repair and exchange when there is a necessary situation of replacing the parts.

### 3. Convenient Filter replacement System

Simple replacement of filters for wavelength changes, no necessity for hand-piece changes reduces trouble when replacement procedure frequently causes mechanical problems.

### 4. Predetermined & Customized programs

- Predetermined programs provide the optimal parameters for various skin problems and skin colors.
- Precisely operated energy and pulse-time assure the superior results.
- Customization of programs is possible for experts.

### 5. Square Pulse Technology system

Due to Square Pulse Technology system, MIPL guarantees stable and duration of energy emission for safe and effective treatment.

### 6. Fast treatment time

MIPL offers fast treatment with no delay and large spot size.

### 7. Other conveniences

- Optimal parameter programs based on a clinical study
- User-friendly 10.4 LCD / touch screen
- Optional hand-piece
- Noiseless

## Innovation User-Interface

### Easy to Use - Preset Program

- Built-in present parameters are optimal for each treatment area to be operated conveniently
- 4 modes/ 4 programs in each mode/ total 16 programs

### Trim and Accurate User Interface

- With easy to use GUI(Graphic User Interface) operators can recognize all buttons intuitively without fail
- LG, 10.4 Inch sensitive touch screen



Intro



Main



Preset



Administrator